promed feeling

Pflegegerät zur Maniküre und Pediküre
Manicure and pedicure care unit
Appareil de soins de manucure et de pédicure
Strumento per manicure e pedicura
Aparato de manicura-pedicura
Verzorgingsapparaat voor de manicure en pedicure
Прибор для ухода за ногтями рук и ног, для маникюра и педикюра
Urządzenie pielęgnacyjne do zabiegów manicure i pedicura

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
FI

GETTING STARTED

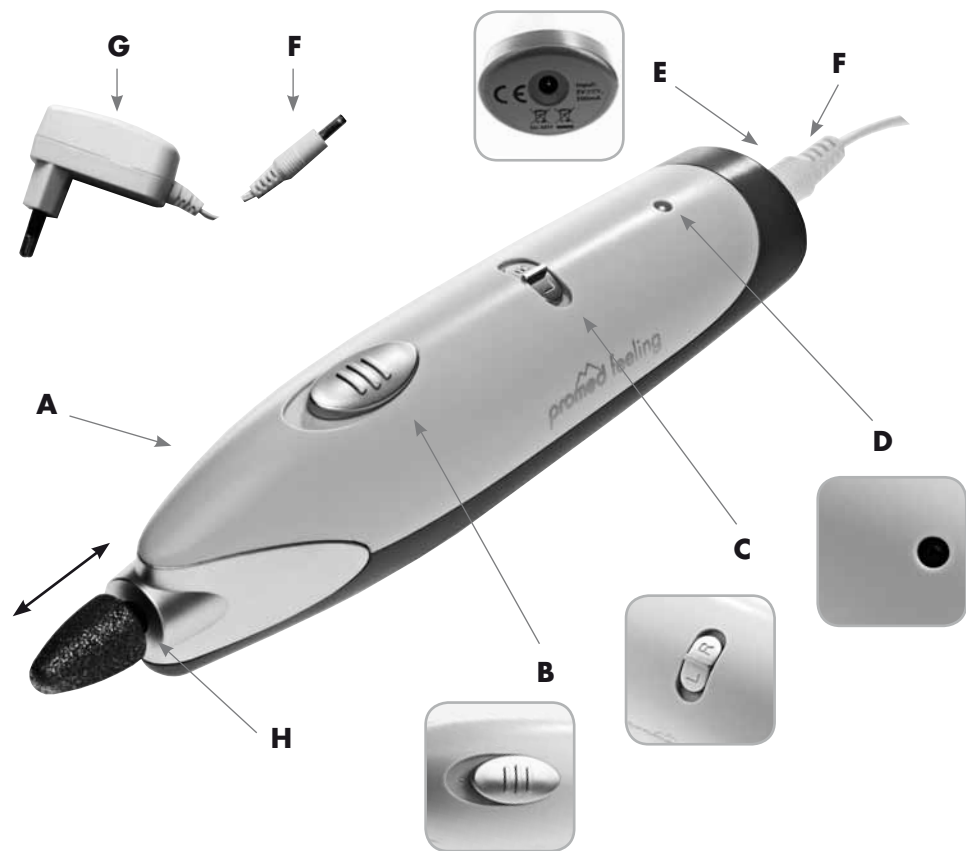- A** Appliance
- B** On/Off Switch, Speed regulator (speed I and II)
- C** Direction switch (L - O - R)
- D** LED-Display
- E** Adaptor socket
- F** Adaptor connector plug
- G** Adaptor
- H** Quick collet

DE

EN

FR

IT

ES

NL

RU

PL

SE

FI

Please read the user manual thoroughly before using the appliance for the first time. Thank you for purchasing the **Promed feeling** manicure and pedicure appliance. You have opted for a German quality product of a leading manufacturer of skin and nail care appliances. Expert advice is available at all places where you can buy **Promed** products or feel free to contact us with your queries. We can recommend a consultant. The equipment was developed specially for the use as manicure and pedicure appliances in the domestic private use.

When you purchase the **Promed feeling** you get a top-quality and elegant manicure/pedicure set, with luxurious ergonomics and professional-quality filing accessories. It is easy to use, making it ideal for hand, foot and nail care, especially for ingrowing, hardened and thickened nails and for treating corns, calluses and hard skin. It is also suitable for cosmetic care for your hands, fingernails and artificial nails. The promed emotion set offers gentle and reliable hand and nail care and is extremely practical to use. The sapphire filing bits are non-wearing and the unit features a safe, automatically-locking chuck designed to be used with a wide range of filing bits and accessories of the standard diameter (\varnothing 2,332 - 2,350 mm).

We, as the manufacturer, cannot be made liable in any way for injury or damages to people or objects that arise from failure to comply with this warning.

We wish you a lot of enjoyment with your new **Promed feeling** unit. We would like to familiarize you with the unit in the following sections. Please read the user manual thoroughly before using the unit for the first time.

Scope of delivery:

Device incl. cable, transformer plug and 5-part manicure. Additional attachments are also available from specialist shops. Inform yourself regarding the details.

INFORMATION ON THE UNIT / ACCESSORIES

Before operating the unit for the first time, please check that all contents have been delivered and that the line voltage is correct. The supplied power adapter (**G**) is an electronic power supply unit intended for an operating voltage of 100 V~ to 240 V~.

- The **Promed feeling** unit uses a high performance motor. The electronic speed control is used to smoothly set the speed of the grinding head. The unit can rotate clockwise or counterclockwise (reversible).
- The grinding head holder is able to fixate all commercial grinding heads with standardized dental size shafts of 2.332 to 2.350 mm (0.091 - 0.093 inches). Please use only high-quality grinding heads with a rounded shaft end and a maximal shaft length of 36 mm (1.26 inches) to guarantee the longest life possible for the grinding head holder!
- Always make sure that the grinding head does not come in contact with other objects when you lay the handpiece down. There is a risk of the shaft becoming compressed and asymmetrical. Asymmetry in the grinding head causes damage to the handpiece!
- Grinding head bits are available at your local specialist retailer or **Promed** (www.promed.de).

DE

EN

FR

IT

ES

NL

RU

PL

SE

FI

INFORMATION ON THE UNIT / ACCESSORIES

- To ensure hygienic operation, the grinding heads must be disinfected after use.
- The handpiece may only be disinfected using disinfectant wipes.
- Information about suitable disinfectants and products can be obtained from **Promed**.

For disinfection of instruments and grinding heads:

Promed Pure ID (Product No. 330810)

For disinfection of surfaces and devices:

Promed Pure FD (Product No. 330815)

For disinfection and cleaning of instruments and grinding heads:

Promed UC-50 (Product No. 330210)

OPERATING THE UNIT

- A.** The appliance functions only without power connection. The batteries must be charged before using the appliance for the first time.

Charging batteries:

Do not recharge the batteries after each use. If the batteries are not used for longer period, it should be first recharged, whenever required. Plug the adaptor connector **(F)** into the adaptor socket **(E)** on the appliance end and connect the adaptor **(G)** to the mains supply. Charging of the batteries takes approx. 16 hours. The LED-display **(D)** lights up to indicate that the charging is in progress. On completion of charging, disconnect the charger from the mains. Once the appliance starts turning slowly, this indicates that the battery must be charged.

- B.** Insert the required abrasive wheel into the relevant quick collet **(H)** of the appliance. The appliance is now ready for operation.

C. Right/Left run

The direction switch **(C)** for the right/left run is located above the LED **(D)**. The corresponding running direction is selected by sliding this switch. The unit is in standby mode in the centre position.

Note:

The direction of run must never be changed during the operation, as the appliance can otherwise be damaged.

D. Speed Regulation

Move the on/off switch **(B)** from "off-(O)" position to I position, to start the manicure. You can further increase the speed by pushing the switch in the same direction at stage II.

- E.** Before using the device for the first time, allow the hand piece to run for approx. 30 minutes with approx. 7,000 rpm (switch position I). This is necessary so that the motor brushes can be grounded.

Note:

A different running noise between clockwise and counterclockwise rotation is normal and caused by bedding in of the carbon brushes or by running in of bearings.

DE

EN

FR

IT

ES

NL

RU

PL

SE

FI

OPERATING THE UNIT

In cases of malfunctioning:

- Check that the power supply is correct.
- Check whether the LED (D) lights red
- Please check whether the rotation direction switch is in the centre position (Standby).
- If the problems have still not been solved after following these steps, please consult your service center immediately.

SPECIAL QUALITY FEATURES

Quick look chuck:

The quick look chuck in the grinding head holder allows you to change grinding heads, easily, quickly and smoothly. The shaft locks automatically as soon as the grinding head starts to rotate.

USAGE

Skin and nail care is not merely a matter that is indicated when suffering correspondingly, it is also a cosmetic matter.

Start with cutting and forming toenails. They absolutely must be filed straight in order to hinder nails growing into the skin. Micro-fine filing using the sapphire disc protects against ingrown nails. The micro-fine files with the sapphire disc prevent an ingrowth. Troublesome calluses always develop in the areas where the foot is over-stressed, mostly the heel, the front balls and the side of the big toe. Treat the callus using the sapphire cone or even better with the drainage grinder. Repeat this procedure regularly. Under no circumstances should you take radical measures when removing the callus. Ingrown toenails can – as needed – be loosened gently using the small cutter. Thick lignified nails can be ground off slowly from the top and filed thinner from the sides.

Crumbly nail areas can also be cleaned using the fine sapphire cutter. Corns are the painful result of a slow callus process. They are calluses with an externally hardened part which is soft inside. Remove the upper callus with the sapphire cone and then work further using the small sapphire cutter. This already removes the pressure. Deeply sunken corns should not be removed by oneself; they should be removed by a doctor. If skin and nail care are carried out systematically using the **Promed** devices and accessories, regular care becomes a positive experience. The special abrasive wheels are described more closely on the following pages. The abrasive wheels are available in various forms and degrees of hardness and can be reordered at any time. Ask your specialist dealer.

Note:

Do not wash your hands or, when applicable, your feet before treatment. A soothing and calming moisturizer should be first used after treatment.

DE

EN

FR

IT

ES

NL

RU

PL

SE

FI

APPLICATION OF THE ABRASIVE WHEELS

Stone sander, fine, pointed

For the removal of difficult corners and edges at those sites that thicker grinding heads cannot access. (Be careful when using the needle burr!)

Stone sander, cylindrical

For the grinding of tough, thickened nails. Thickened nails can be shaped easily and their appearance improved. Pressure is relieved. It is also suitable for grinding the edges of hard corns to relieve pressure. (Only grind away the areas of pressure!)

Sapphire cutter, long

To remove hard corners and edges. For the treatment of nail underside, for example removing dirt.

Sapphire file, cylindrical

For the grinding of tough, thickened nails. Thickened nails can be shaped easily and their appearance improved. Pressure is relieved. It is also suitable for grinding the edges of hard corns to relieve pressure. (Only grind away the areas of pressure!)

Grinding cone, sapphire

For callused and rough skin on the feet and hands, but also for abrading very thick, hard nails. It reshapes crooked nails and relieves painful pressure.

DE

EN

FR

IT

ES

NL

RU

PL

SE

FI

MAINTENANCE AND CARE

Cleaning the unit

- Regular maintenance and cleaning ensure optimal results and extends the life of your unit.
- Store the **Promed feeling** in a safe, cool and dry location when it is not in use.
- To ensure that the safety, reliability and performance of the unit are not adversely affected, enhancements, readjustments, alterations or repairs are not permitted aside from those performed by **Promed** authorized persons.
- Before cleaning the unit, make sure that it is turned off and disconnected from the electricity.
- Use a slightly moist, soft cloth to clean the unit.
- Do not use any aggressive substances such as scouring agents or other detergents as they may corrode the surface.
- Wipe clean with a soft, moist cloth and do not use scouring agents.
- Do not immerse the unit in fluid and keep the unit away from solvents and strong detergents.
- Do not repair the unit. The unit does not have parts that consumers can service.

Cleaning the grinding heads

The abrasives can be cleaned from dust and skin residuals by using a cleaning brush (e.g. **Promed cleaning brush** Art.No. 198910). It is always advisable to disinfect the abrasives with a commercially available instrument disinfection (e.g. **Promed PURE-ID**), in order to avoid transmission of infections. Please make sure to keep the pins of the abrasives free of grease and clean these from time to time.

Caring for the rechargeable batteries / precautionary measures when charging

- To maintain full battery capacity, it is recommended that only fully discharged batteries be charged again. Disconnect the power supply unit from the electricity network when not in use for longer periods of time.
- Do not continuously charge the unit for more than 3 days.
- Do not operate the unit with the power supply unit when the batteries are fully charged.
- Every few months, operate the unit until the motor stops altogether. Then place the switch onto "off" and fully charge the unit up again. Do not charge the unit or store it when the room temperature is under 41° F (5° C) or above 122° F (50° C).
- Do not charge the unit or store it near heaters or other sources of heat or in places where the unit is subject to strong sunlight.
- Do not charge it in moist places.
- Pull the plug of the power supply unit out of the wall socket when you have finished charging.

Changing the battery (service)

Should you need to replace the battery, please send your **Promed feeling** unit to either the manufacturer (**Promed GmbH**) or an authorized Promed customer service center.

DE

EN

FR

IT

ES

NL

RU

PL

SE

FI

SAFETY INSTRUCTIONS

Please pay attention to the following before first time operation:

Never place the unit in direct contact with fire, gas, oxygen or hot objects such as hot plates.

Please observe the following during operation:

A. Grinding heads to be used:

- Standardized dental size: 2.332 - 2.350 mm / 0.091 - 0.093 inches.
- Only use centered grinding heads with a rounded shaft end.
- The use of rusted or bent grinding heads is not permitted.

B. Switch the unit off before you change the direction of rotation.

C. Appliance power supply: 100 V~ till 240 V~.

D. Switch the unit off before you change the grinding head.

E. Take every precaution possible to ensure that the unit does not fall down or get damaged in any other way.

F. If problems do arise with the handpiece, take it immediately to get it repaired.

G. Do not lubricate or wash the handpiece.

H. For home use only.

I. Scrupulous supervision is required when the device is being used in the vicinity of children.

J. Protect the battery against heat, e.g. also against continuous sun exposure and fire. It can cause explosion.

K. If you do not use or charge the appliance, disconnect the adaptor from power supply.

DANGER!

- Do not place the unit in direct contact with water or other fluids.
-  Never place or use the unit in a wet or moist environment.
- Do not touch the unit with wet hands.
- Operate the appliance only with an appropriate power adapter.
- It is not permitted to use the unit outdoors.

WARNING!

- Make sure that children do not play with the device.
- Never leave the unit unattended when children or inexperienced persons are in its vicinity.
- Only use the unit for the applications described in this manual.
- Only use attachments and accessories from the manufacturer.
- Never use the unit after the electrical lead has been damaged, when the unit has fallen down or when it has come in contact with water.
- This unit is not designed to be used by people (including children) who have limited physical, sensory, or intellectual abilities or who do not have sufficient experience and/or sufficient knowledge. This does not apply if these people are supervised by a person responsible for their safety or receive instructions from them on how to use this unit.
- Never let children play with the packaging material; there is a risk of suffocation.
- An electrical appliance is not a toy. Children do not recognize the risks that could develop when using electrical appliances. Therefore, operate and store the unit out of the reach of children.
- Do not place or use the unit on hot surfaces or near naked gas flames to ensure that the casing does not melt.

DE

EN

FR

IT

ES

NL

RU


PL

SE

FI

SAFETY INSTRUCTIONS

- Never let the unit come into contact with water. Use a moist cloth to clean the unit. Do not use solvents.
- Never carry or pull the unit by the power cord. Use the plug to pull the power cord out of the wall socket. Do not pull on the power cord itself.
- Do not pull the power cord over sharp edges, wedge it in between objects, hang it down or position or hang it over hot objects (e.g., hot plates) or naked flames. Protect the power cord from heat or oil.
- Only use extension cords when you have satisfied yourself beforehand that they are not defect or damaged.
- If the power cable of the appliance gets damaged, have it replaced by an expert immediately
- Do not store the unit near a wash basin or bath tub as there is the danger that the unit may fall or be pulled into the wash basin or the bath tub.

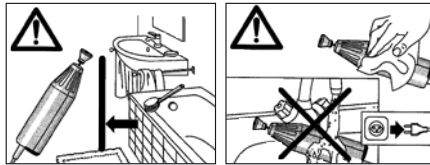- To avoid all risks of injury, blockage or damage to the unit, make sure that the unit does not come into contact with clothing, cords, shoelaces, bristles, wires, threads, hair etc. when in operation.

Note:

The AC adapter may become slightly warm. This does not indicate a malfunction. The device can be used easily.

DE

EN

FR

IT

ES

NL

RU

PL

SE

FI

DISPOSAL

Power tools, accessories, and packaging should be recycled in an environmentally suitable manner.

Only for EU countries:

Do not dispose of power tools with the domestic waste! In accordance with the European Directive 2002/96/EC on waste electrical and electronic equipment and its implementation into national law, nonfunctioning power tools must be collected separately and recycled in an environmentally suitable manner.

This symbol indicates that disposal of this product in the domestic waste is not permitted within the EU. Waste devices contain valuable recyclable material that should be recycled. Recycling also prevents uncontrolled waste disposal from damaging the environment and human health. Therefore, please dispose of all waste devices in the appropriate collection systems or send the device to the place where you bought the device for disposal. They will then recycle the device.

Rechargeable batteries/batteries:

Do not throw your rechargeable batteries/batteries in the domestic waste, in fire, or in water. Rechargeable batteries/batteries should be collected, recycled, or disposed of in an environmentally friendly manner.

Only for EU countries:

According to the Directive No. 91/157/EEC, defect or used rechargeable batteries/batteries must be recycled. Waste rechargeable batteries/batteries can be directly handed to: **Promed GmbH**, Lindenweg 11, D-82490 Farchant.

WARRANTY

This product was produced with the utmost care and was thoroughly tested before it left the factory. Therefore, we provide a warranty period of 24 months from the date of purchase for this product. Promed products conform to their description and the relevant specifications; it is your responsibility to ensure that the product which you have purchased is designed for your intended use.

- In the event of material defects and/or manufacturing errors which appear during proper use and are detected during the warranty period, we will replace all defective parts of the product free of charge within the warranty period, including the labour costs for warranty repairs.
- The warranty does not cover:
 - Normal wear on the product
 - Defects which were caused during transportation or storage of the product
 - Defects or damage which were caused through improper use or inadequate maintenance
 - Damage resulting from disregarding the notes in the operating manual
 - Damage resulting from modifications to the product which were not carried out by Promed.
 - Damage caused by sharp objects or as a result of torsion, compression, being dropped, an abnormal impact or other actions which are beyond Promed's reasonable control.
- Parts which are subject to wear (e.g. moving parts such as ball bearings, etc., fasteners, etc.) are generally excluded from the warranty.
- A claim under the warranty is not enforceable if:

WARRANTY

- the product is not sent back in its original packaging or in adequately secure packaging;
- it has been modified or repaired by a person or company other than Promed or one of Promed's authorised dealers;
- the product has been repaired using parts which are not approved by Promed;
- the product's serial number has been removed, erased, modified or made unreadable;
- For reasons of hygiene, any products which have been exposed to direct contact with a body or with bodily fluids (e.g. blood) must be packaged in an additional plastic bag before they are sent back. In this case, a specific note of this situation must be made in the package or in the accompanying documents.
- For components or products which are replaced during the warranty period, a warranty is only provided for the remainder of the original warranty period; provided that the product or part replacement was carried out by Promed or one of Promed's authorised dealers.
- The warranty period begins on the day of purchase. Claims under the warranty must be established within the warranty period. Complaints arising after the expiry of the warranty period cannot be considered.
- Under these warranty conditions, the warranty shall only apply if the date of purchase can be proven by means of a sales receipt or similar.
- We reserve the right to make technical and visual changes as well as changes to the features.
- This warranty is only valid and enforceable in the country in which it was purchased by the original buyer, provided that it was Promed's intention that the product should be offered for sale in this country. This warranty is also enforceable in any country within the European Economic Area in which Promed has an authorised importer or sales partner. Depending on each country, specific and differing warranties and guarantees may be applicable owing to the particular applicable legislation. These legal provisions are neither excluded nor limited by these warranty conditions. To the extent permissible under national law, the warranty period shall not be extended, renewed or otherwise affected by any subsequent resale, repair or replacement of the product.
- The provisions of the UN Convention on Contracts for the International Sale of Goods shall not apply.
- Our warranty conditions do not affect the seller's legal warranty obligation.
- To the fullest extent possible under the applicable mandatory legislation, these warranty conditions constitute your sole and exclusive legal remedy and shall apply in place of any other express or implied warranty conditions. Promed is not liable for any unusual, incidental, punitive or consequential damages including, but not limited to loss of profits, loss of use, loss of revenue, costs of replacement equipment or facilities, third party insurance claims, and damage to property which are the result of the purchase or use of the product or occur owing to a breach of warranty, a breach of contract, negligence, product defects or any other legal or statutory circumstances, even though Promed was aware of the possibility of such damages. Promed is not liable for any delay in warranty claims.
- Promed cannot be held liable for any possible translation errors.

The following information is essential for smooth processing:

1. Original proof of purchase/receipt or dealer's stamp with date of purchase
2. Defect discovered
3. Product designation / type / serial/batch number

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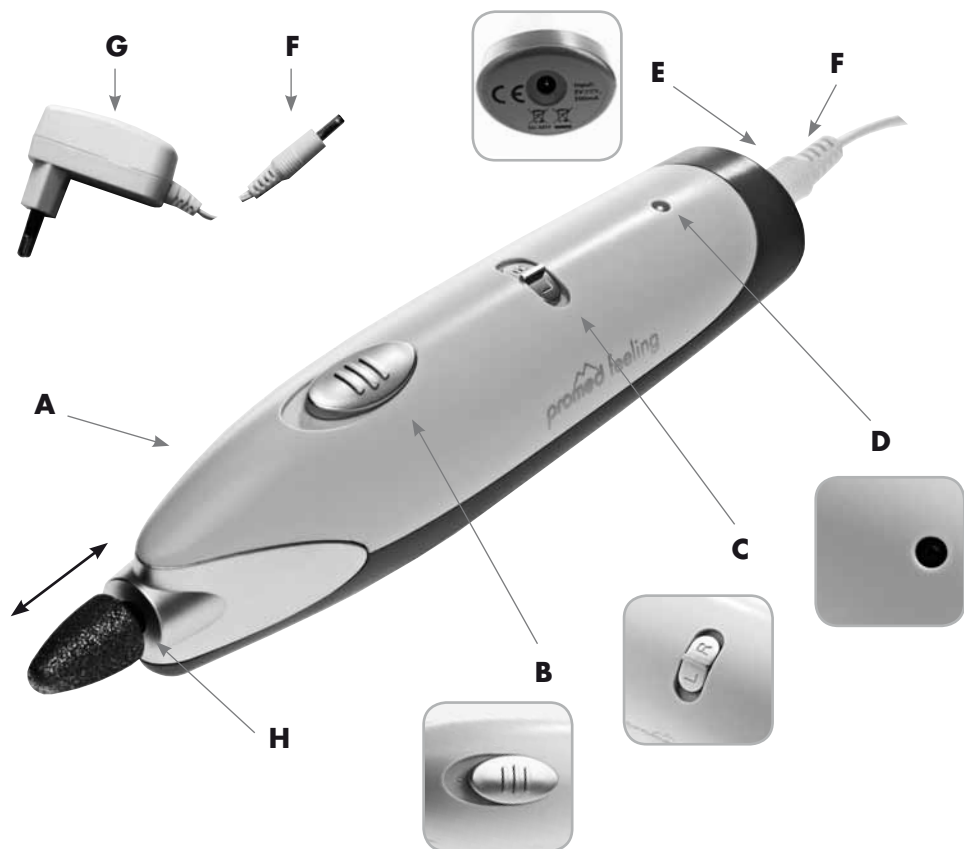
Symbol for protection class II appliances

The symbol corresponds to symbol 5172 of the EN 60417-1:2003 Standard. The appliance is double insulated as a protective measure. This prevents conductive parts of the appliance from becoming impaired, something that could lead to malfunctioning of the basic insulation voltage.

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Symbol for WARNING! Indicates sections of the user manual that are relevant for safety!

The symbol corresponds to symbol 0434 of the DIN ISO 7000:2008 Standard. Instructions for an appliance, system or accessory that contain important information for the user that are relevant for safety.

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Pay attention to the symbol for the user manual!

The symbol corresponds to the symbol 1641 in the DIN ISO 7000:2008 Standard. It is imperative that the manual be read before using the unit.

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EN | The power supply unit complies with the Ecodesign Directive (ErP) No. 2013/617/EC requirements for external power supply units and the requirements of phase 1 and 2 of the EC Directive 2013/617/EC.

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The logo for 'promed feeling' features a stylized mountain range icon above the brand name. 'promed' is in a bold, lowercase sans-serif font, and 'feeling' is in a lighter, lowercase sans-serif font.

Pflegegerät zur Maniküre · Care device for manicure · Appareil de soin pour manucure · Apparecchio per la manicure · Aparato de manicura · Verzorgingstoestel voor manicure · Маникюрный прибор · Urządzenie pielęgnacyjne do manicure · Manikürmaskin · Manikyyri-laite

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